URIT-8210

Automatic Chemistry Analyzer

Specifications

Assay methods: End Point, Kinetic, Fix time etc.

Principle: Photoelectric colorimetry
Light Source: Halogen lamp 12V/20W

Photometry range: 0~3.2Abs Resolution: 0.0001 Abs

Wavelength: 10 wavelengths optional (340nm~800nm)

Throughput: 330tests/hour

Reagent tray: 59 reagent positions, 1 detergent position

Sample tray: 71 sample positions, including detergent, standard, QC, STAT positions

Reaction tray: 90 reaction cuvettes

Sample volume: 2~100ul, with 0.1ul increment

Reagent volume: R1:10~500uL, R2: 10~500uL, with 0.5uL increment

Minimum reaction volume: 150uL Maximum reaction time: 10 minutes

Water Consumption: 6L/hour under working status

Clean unit: 8-step auto-washing system with detergent and 37℃ water

Calibration: Calibration reset, select best test point by reaction curve, no need second calibration

Line/non-line;multi-standards assay

Control rules: Westgard multi-rule, Cumulative sum check, Twin plot

3 level controls for each ifem, analyxing and printing QC analysis diagram

Temperature control: incubator 37 ± 0.1℃

Power supply: AC110/230(1 ± 10%)V,50/60Hz,500VA

Ambient: Operating temperature: 10~30℃

Relative humidity: ≤85%

Atmospheric pressure: 86~106kPa

Features

- ▼24 hour non-stop cooling system to ensure reagent at 2-8℃
- ▼ Dual-unattached durable ceramic syringes to ensure accuracy & precision
- Collision protection in both vertical and horizontal directions, stop & alarm automatically once touching barrier, not affect former tests
- Automatic eligible cuvettes detection & selection
- Select best test point by reaction curve, create new factor automatically
- High accurate optical system
- ▼ Support LIS interface

URIT MEDICAL ELECTRONIC (GROUP) CO., LTD.


Head Office Add: No.4 East Alley, Juhus Road, Gullin, Guangxi 54 100 1, P.R. China Tat: +66-773-2206566 Fax: +66-773-2206560

Shenzhen Office Add.: Room EJ, 24h Floor, New Bachul Building, Nanhai Road Nanshan Datrict, Sherothan, P.R. China Tat. +68-755-20030833 Fax: +68-755-86038150





Global Diagnostics Supplier

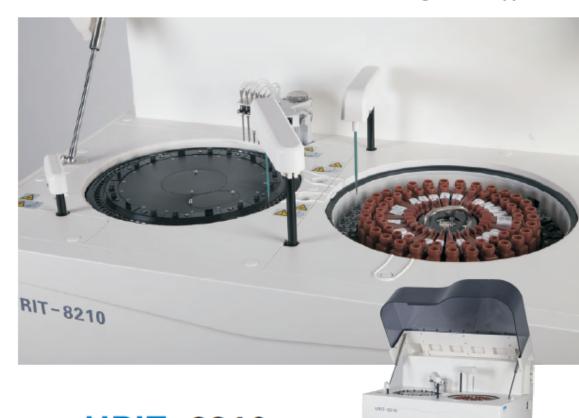

URIT-8210

Automatic Chemistry Analyzer

- Random access, 330 tests/hour
- ISE and barcode (Optional)
- 10-wavelength, high accurate and closed optical system
- Unattached probe and mixer
- Supports LIS system

